

ACROSS CTO™ OTW


by Acrostak

The world's smallest lesion entry profile

Designed to cross the most challenging CTOs thanks to superior shaft technology.

Outstanding Balloon performances

- The world's smallest lesion entry profile of 0.015" (0.38mm)
- Truly tapered tip with less than 2° inclination
- Unsurpassed material resistance with **RBP of 17 atm**
- **Hydrolubric™ coating** for minimal friction and optimal trackability

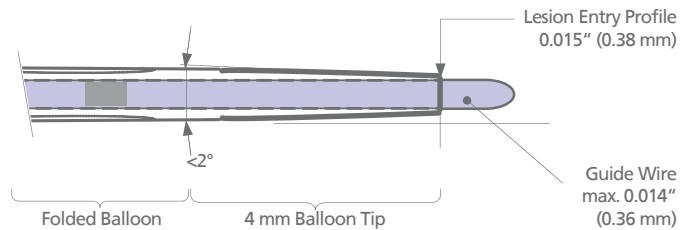


Hypotube shaft design


- Metallic hypotube with ultra thin shaft wall and extremely smooth surface
 - ▶ **Provides extra push** to penetrate very stiff occlusion caps

Flexy shaft design


- **Proximal shaft made out of flexible polymer material**
 - ▶ **Provides extra flexibility** to navigate within very tortuous anatomies, often encountered in **retrograde approach**

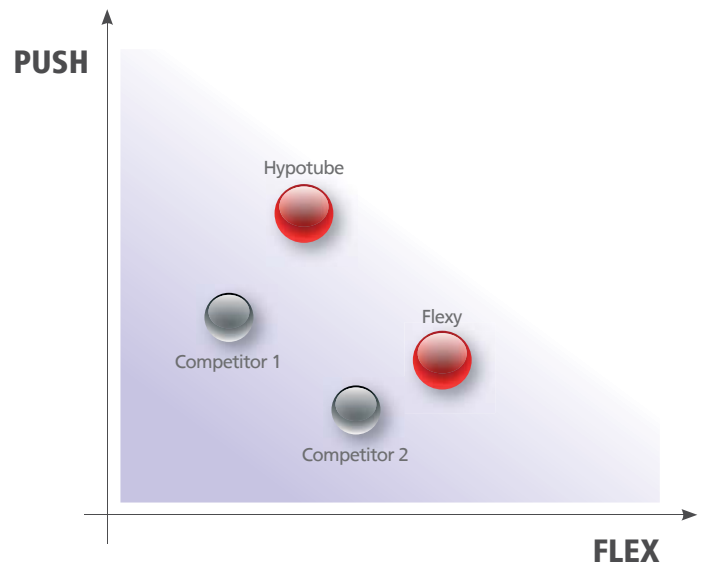

ACROSS CTO™ OTW (hypotube)

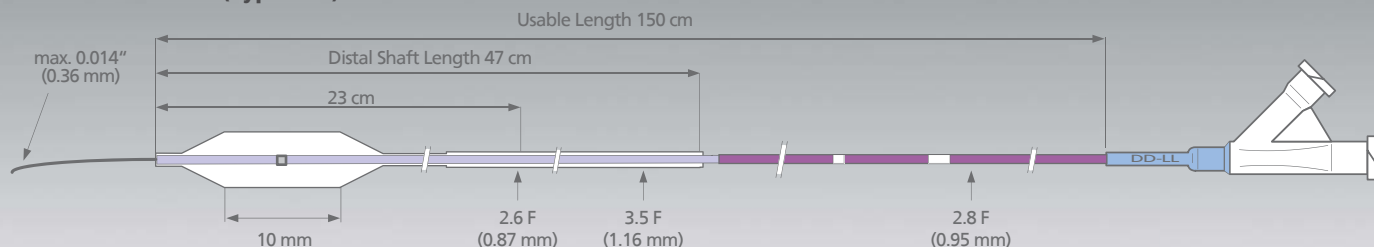

ACROSS CTO™ OTW flexy

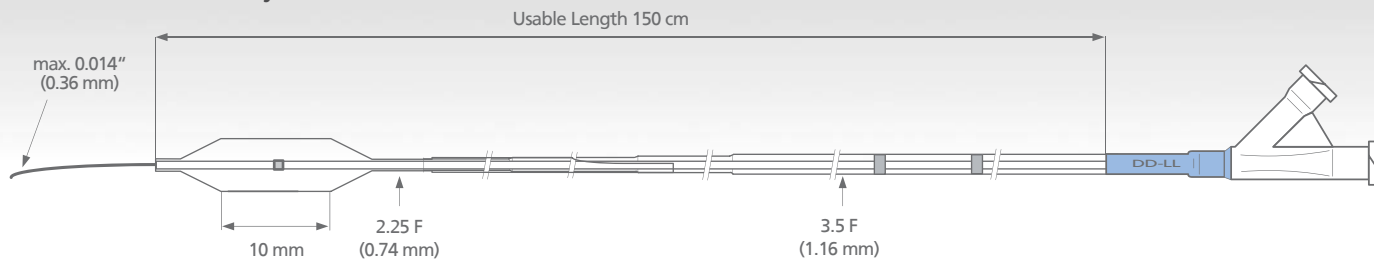

COMPLIANCE CHART

Pressure (atm)	Ø1.1 10/15/20
5	1.04
7 (709 kPa)	1.10 ^N
9	1.13
11 (1115 kPa)	1.16
13	1.20
15	1.23
17 (1722 kPa)	1.28*
19	1.38
21	1.45

^N Nominal Pressure * Rated Burst Pressure. Do not exceed.

ORDERING INFORMATION

Reference no.	Diameter (mm)	Length (mm)
ACROSS CTO™ OTW*		
110100481	1.1	10
110150481	1.1	15
110200481	1.1	20
ACROSS CTO™ OTW flexy		
110100482	1.1	10
110150482	1.1	15
110200482	1.1	20

* Hypotube



Worldwide distribution by:

Acrostak International
Route de Malagnou 6
CH-1206 Geneva
SWITZERLAND

sales@acrostak.com
Phone +41 22 347 68 70
Fax +41 22 347 68 71

Manufactured by:

Acrostak (Schweiz) AG
Stegackerstrasse 14
CH-8409 Winterthur
SWITZERLAND

info@acrostak.com
Phone +41 52 233 95 51
Fax +41 52 233 95 50

www.acrostak.com

Distributed by:

