


GRIP™
by Acrostak

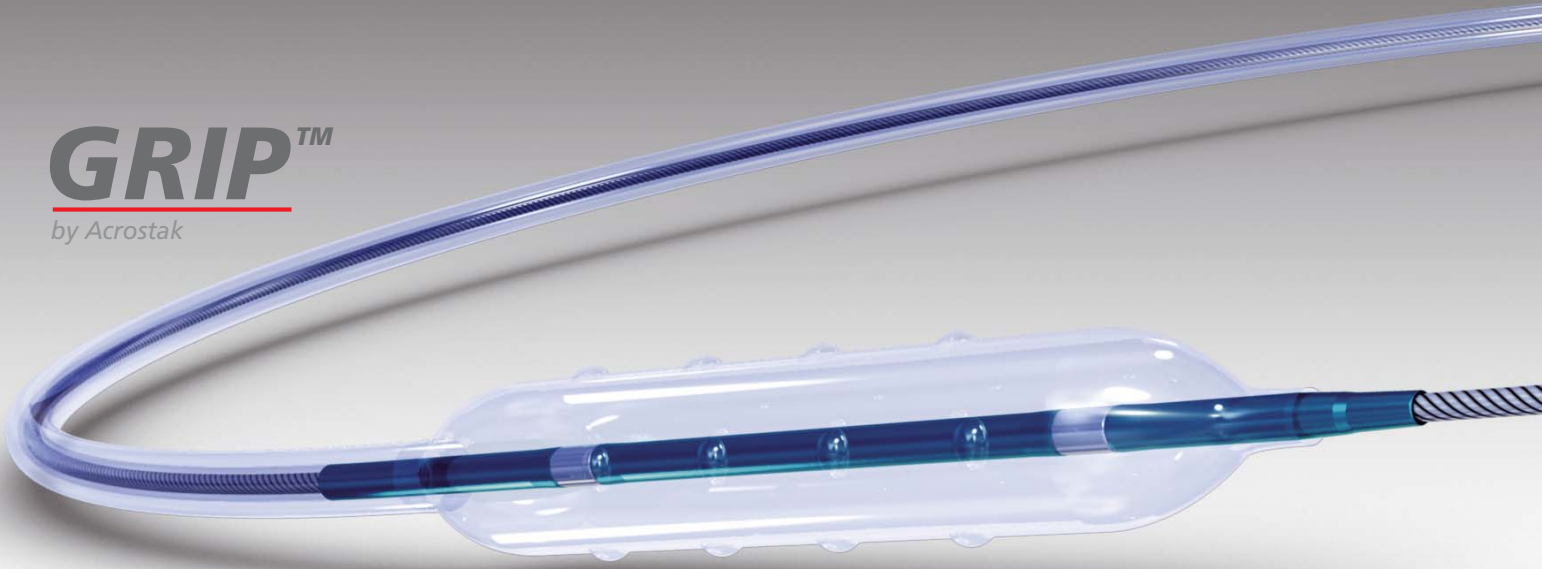

Full control over slippery and calcified lesions


Optimal for the treatment of **intra-stent restenosis, heavily calcified lesions** and for **stent post-dilation**.

Unique Knobs Design

- 4 lines of knobs on the surface of the balloon anchor the balloon to the lesion during balloon inflation
 - ▶ **Reduced risk of dissection** caused by balloon slippage at the stent margins of ISR, and in delicate areas such as bifurcations, ostial lesions, and short calcified lesions
- 4 lines of knobs designed to crush the hard plaque remaining between the struts of the stent
 - ▶ **Optimal plaque crushing** for superior post-dilation

Outstanding Material Strength

- High pressure balloon with **RBP of 22 atm** for all sizes
- **Ultra-resistant to puncture**
- **Truly tapered tip:**
short tip design for maximum push (GRIP™)
long tip design for maximum trackability in tortuous anatomies (GRIP™TT)
- ▶ **Securely reach and dilate** heavily calcified lesions without compromising **crossability**

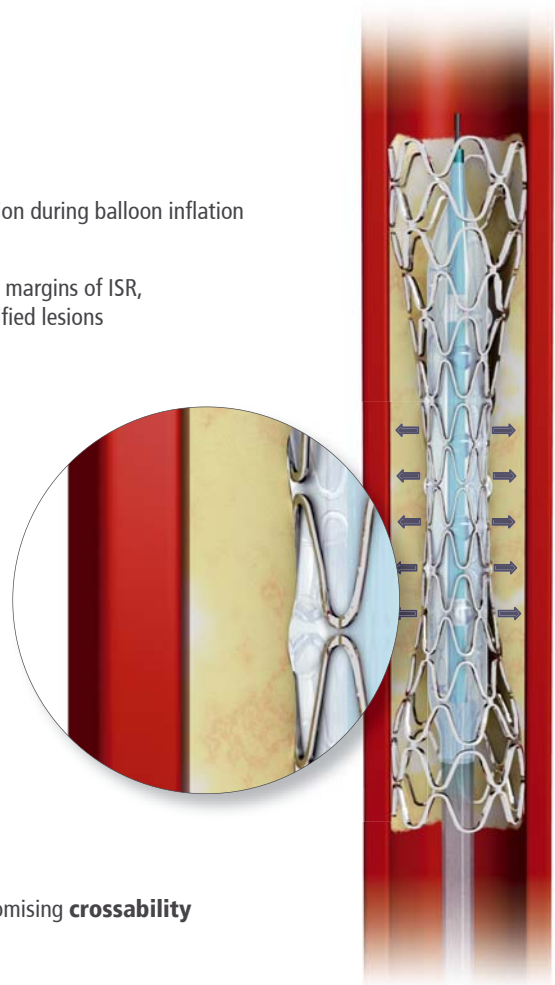


CE 0481


© ACROSTAK | All rights reserved.

Acrostak
Devices for interventional applications

GRIP™

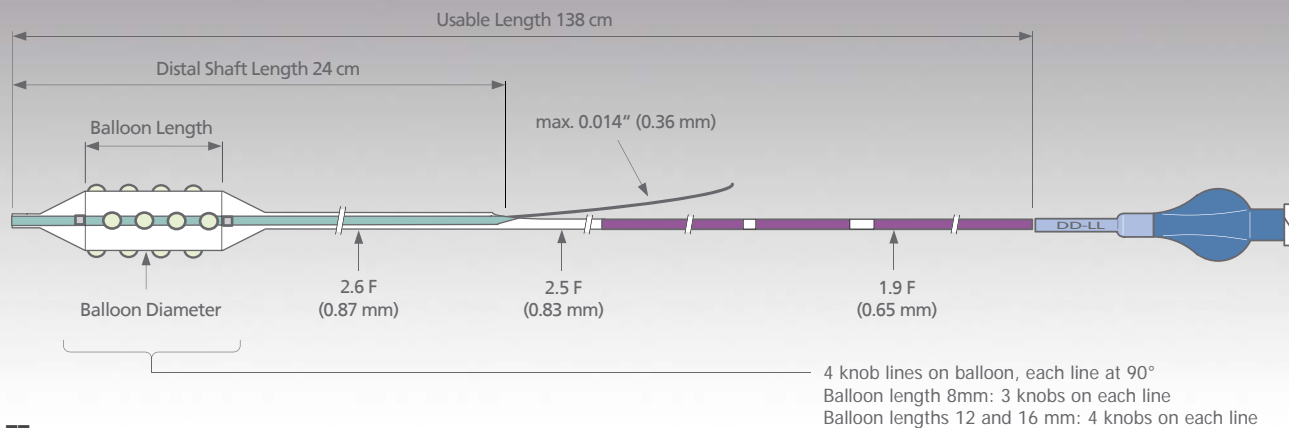

GRIP™ TT

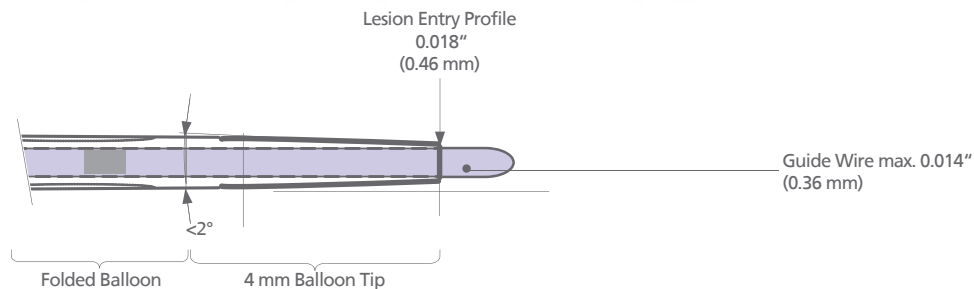

COMPLIANCE CHART (GRIP™ & GRIP™ TT)

Pressure (atm)	Ø2.5	Ø3.0	Ø3.5	Ø4.0
	10/15/20	10/15/20	10/15/20	10/15/20
6	2.30	2.76	3.19	3.65
8	2.37	2.84	3.31	3.78
10	2.45	2.94	3.44	3.91
11 ^N (1115 kPa)	2.50	3.00	3.50	4.00
12	2.54	3.05	3.56	4.07
14	2.62	3.14	3.67	4.19
16	2.68	3.22	3.76	4.29
18	2.73	3.27	3.83	4.37
20	2.77	3.33	3.90	4.45
22 [*] (2229 kPa)	2.81	3.38	3.97	4.52

* Rated Burst Pressure - do not exceed

^N Nominal

ORDERING INFORMATION

Reference no.	Diameter (mm)	Length (mm)
GRIP™	GRIP™ TT	
250080340	250080660	2.5 8
250120340	250120660	2.5 12
250160340	250160660	2.5 16
300080340	300080660	3.0 8
300120340	300120660	3.0 12
300160340	300160660	3.0 16
350080340	350080660	3.5 8
350120340	350120660	3.5 12
350160340	350160660	3.5 16
400080340	400080660	4.0 8
400120340	400120660	4.0 12
400160340	400160660	4.0 16



Worldwide distribution by:

Acrostak International
Route de Malagnou 6
CH-1206 Geneva
SWITZERLAND

sales@acrostak.com
Phone +41 22 347 68 70
Fax +41 22 347 68 71

Manufactured by:

Acrostak (Schweiz) AG
Stegackerstrasse 14
CH-8409 Winterthur
SWITZERLAND

info@acrostak.com
Phone +41 52 233 95 51
Fax +41 52 233 95 50

www.acrostak.com

Distributed by:

